
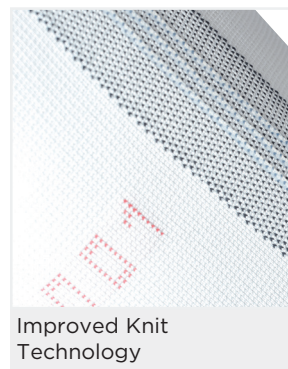


The Streifeneder silicone liners offer solutions to meet a wide range of user requirements. Each liner uses robust, first-class materials, alongside Streifeneder's new knit technology, to provide a great balance of comfort and function.


The PrimoSil.basic is a 3 mm silicone liner made with especially soft silicone.

Features and Benefits

- Ideal silicone liner for those with pressure sensitive residual limbs
- Suitable for low to moderate activity levels
- Designed for those looking for additional comfort and compliance
- No antibacterial finish for users with more sensitive skin
- Integrated 10cm knitted matrix reduces distal elongation under tensile load
- Easy to don and doff due resilient polyamide/elastane cover
- Integrated coloured vertical stripe for easier donning orientation
- Distal silicone-rotation control
- Available with or without distal connections

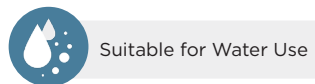

PrimoSil.basic

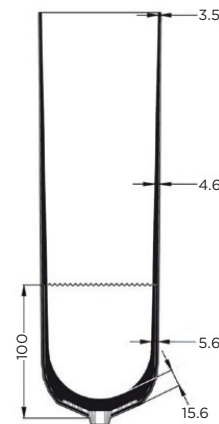

Ordering Information

PrimoSil.basic				
Ordering Example NMA44L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
16, 18, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	10cm	with	NMA44L10/Size
16, 18, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	without	without	NMA44L20/Size

**Sizing Note:** Measure the distal circumference approx. 4 cm above the distal end and select the liner in the same size or one size smaller.

Suitability Chart

Amputation Type	PrimoSil.basic
Pear-shaped below-knee (Transtibial)	● 3mm basic
Normal, cylindrical below-knee (Transtibial)	●●● 3mm basic
Bony, conical below-knee (Transtibial)	●●● 6mm soft

**Key:** Optimal: ●●● Good: ●● Suitable: ●

Technical Information

Activity Level	K1 to K3
Wall Thickness	3mm
Total Length	Approx 370mm
Max. Patient Weight	150kg
Material	Especially soft silicone, medical grade, approx. 26 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter or without
Warranty	6 Month
Colour	Sand

Supplied with a drying stand.